


**An Ultra Aqueous C18 column enhances retention of pyridine while maintaining a high level of base deactivation.**

Ultra Aqueous C18 (cat.# 9178565)

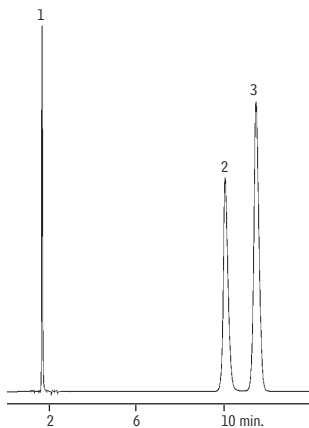

***\*Enhanced selectivity & excellent peak shape!***


**Peak List:**  
1. uracil  
2. pyridine  
3. phenol

**Column:**  
Dimensions: 150 x 4.6mm  
Particle size: 5 $\mu$ m  
Pore size: 100Å

**Conditions:**  
Mobile phase:  
20mM potassium phosphate pH 7.0:  
acetonitrile (80:20)

Flow: 1.0mL/min.  
Temp: 25°C  
Det.: UV @ 254nm

Conventional base-deactivated C18

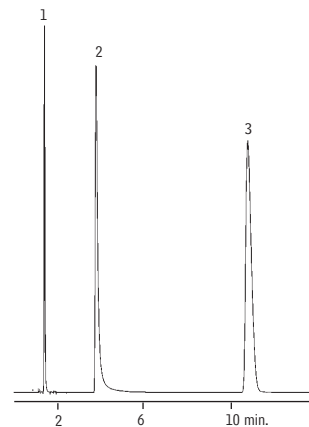

LC\_0143