


# Barbiturate Drug Panel on Raptor™ Biphenyl by LC-MS/MS

Peaks	tr (min)	Precursor Ion	Product Ion	Product Ion
1. Phenobarbital	1.97	230.9	187.8	85.0
2. Butalbital	2.21	223.0	180.0	84.9
3. Amobarbital	2.95	225.0	182.0	84.8
4. Pentobarbital	3.04	225.0	182.0	84.9
5. Secobarbital	3.69	237.0	193.9	84.8

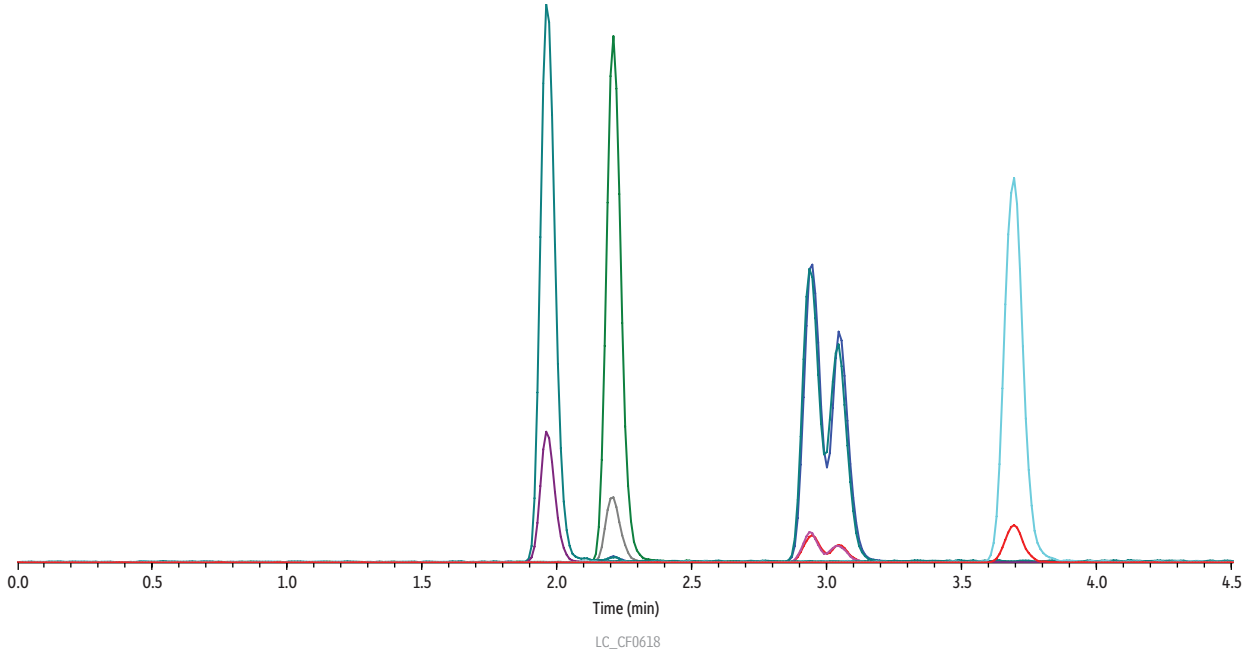

**Column** Raptor™ Biphenyl (cat.# 9309A12)  
 Dimensions: 100 mm x 2.1 mm ID  
 Particle Size: 2.7 µm  
 Guard Column: Raptor™ Biphenyl EXP® guard column cartridge 5 mm, 2.1 mm ID, 2.7 µm (cat.# 9309A0252)  
 Temp.: 30 °C

**Sample**  
 Diluent: Water  
 Conc.: 500 ng/mL  
 Inj. Vol.: 5 µL

**Mobile Phase**  
 A: Water + 0.1% formic acid  
 B: Methanol + 0.1% formic acid

Time (min)	Flow (mL/min)	%A	%B
0.00	0.6	55	45
4.00	0.6	45	55
4.01	0.6	55	45
6.00	0.6	55	45

**Detector** MS/MS  
**Ion Mode:** ESI-  
**Mode:** MRM  
**Instrument** UHPLC