
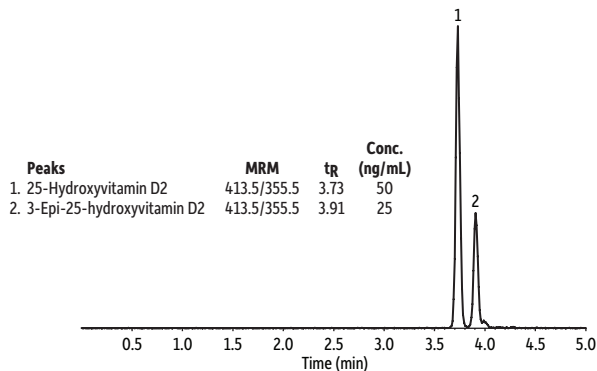


25-Hydroxyvitamin D2/D3 and 3-Epi-25-Hydroxyvitamin D2/D3 on Raptor™ FluoroPhenyl by LC-MS/MS

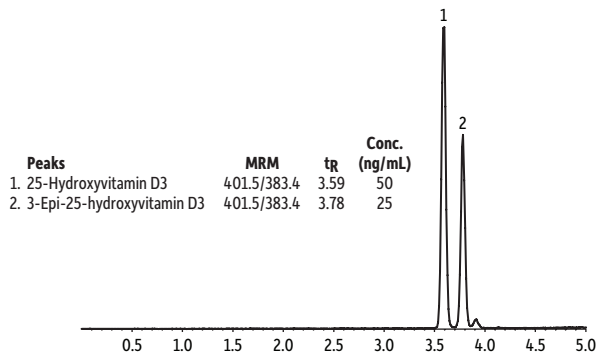

LC_CF0644

Column	Raptor FluoroPhenyl (cat.# 9319A1E)		
Dimensions:	100 mm x 3 mm ID		
Particle Size:	2.7 µm		
Temp.:	30 °C		
Sample			
Diluent:	Water:methanol (50:50)		
Conc.:	25-50 ng/mL		
Inj. Vol.:	5 µL		
Mobile Phase			
A:	0.1% Formic acid in water		
B:	Methanol		
Time (min)	Flow (mL/min)	%A	%B
0.00	0.6	25	75
4.00	0.6	15	85
4.10	0.6	0	100
5.00	0.6	0	100
5.01	0.6	25	75
7.00	0.6	25	75

Detector MS/MS
Ion Mode: ESI+
Mode: MRM
Instrument HPLC