


## SRM 870 Test Mix on Viva C18

1. uracil
2. toluene
3. ethylbenzene
4. quinizarin
5. amitriptyline

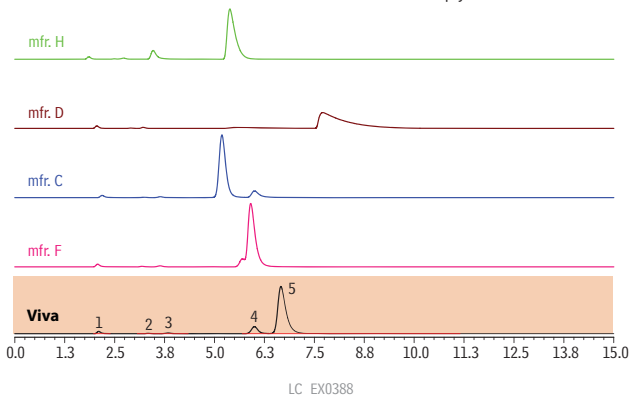

### Sample:

Inj.: 10 $\mu$ L HPLC Performance (NIST 870) Test Mix (cat.# 31699)

### Column: Viva C18

Cat.#: 9514562  
Dimensions: 150 x 2.1mm wide pore  
Particle size: 5 $\mu$ m  
Pore size: 300 $\text{\AA}$

### Conditions:

Mobile phase: water:methanol  
(25:75 v/v)  
Flow: 0.2mL/min.  
Det.: UV @ 254nm