


Lot-to-Lot Reproducibility of Allure Acrylamide

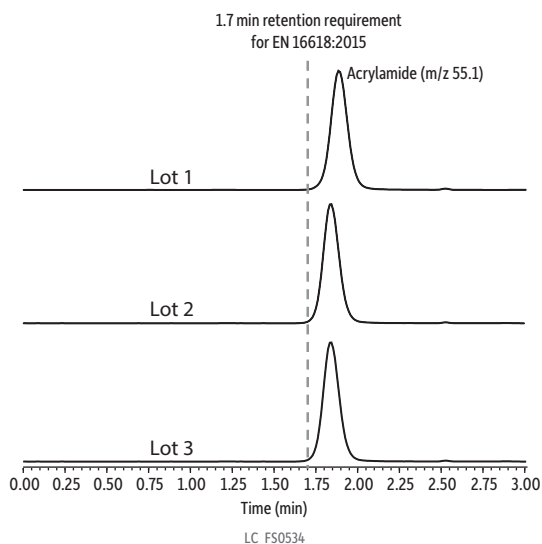

Peaks	(ng/mL)	Conc. Precursor	Product
1. Acrylamide-d3 (IS)	200	75.1	58.1
2. Acrylamide	200	72.1	55.1

Column Allure Acrylamide (cat.# 9167552)
Dimensions: 50 mm x 2.1 mm ID
Particle Size: 5 µm
Pore Size: 60 Å
Guard Column: Allure Acrylamide 10 mm, 2.1 mm ID, 5 µm (cat.# 916750212)
Temp.: 22 °C
Sample Acrylamide (cat.# 30494)
Diluent: Water
Inj. Vol.: 10 µL
Mobile Phase
 A: 0.001% Formic acid in water
 B: 0.001% Formic acid in acetonitrile

Time (min)	Flow (mL/min)	%A	%B
0.00	0.4	100	0
1.00	0.4	100	0
2.00	0.4	10	90
2.01	0.4	100	0
5.50	0.4	100	0

Detector MS/MS
Ion Mode: ESI+
Mode: MRM
Instrument HPLC