


Retention Time Stability for Roc Cyano

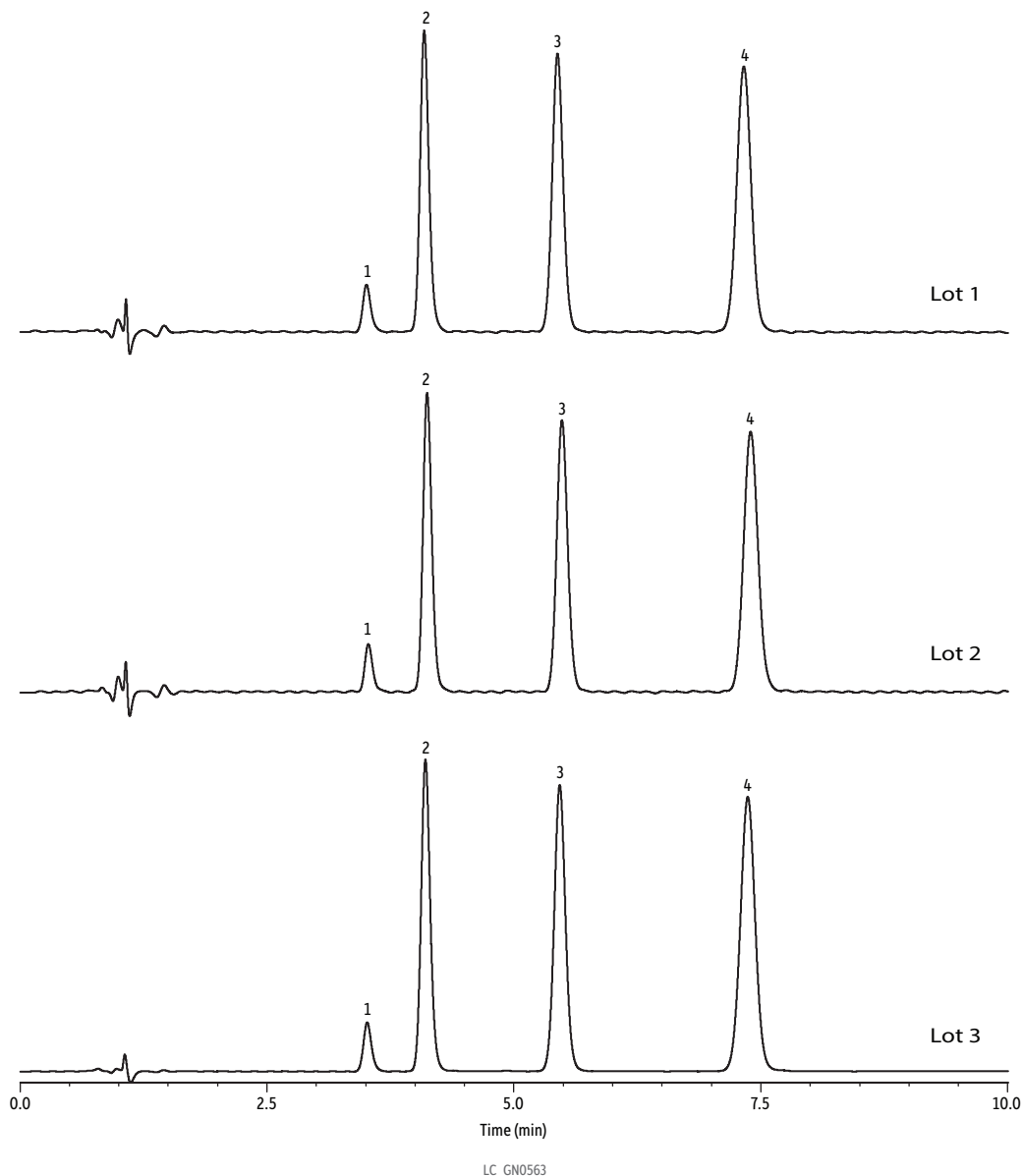

Peaks
1. Benzene
2. Nitrobenzene
3. 1,3,5-Trinitrobenzene
4. 2,4-Dinitrotoluene

tr (min)
3.5
4.1
5.5
7.4

Column	Roc Cyano (cat.# 953631E)
Dimensions:	100 mm x 3.0 mm ID
Particle Size:	3 μ m
Pore Size:	100 \AA
Temp.:	40 $^{\circ}$ C
Inj. Vol.:	1 μ L
Mobile Phase	Water:acetonitrile (80:20)
Flow:	0.6 mL/min
Detector	UV/Vis @ 254 nm
Cell Temp:	40 $^{\circ}$ C
Instrument	HPLC