
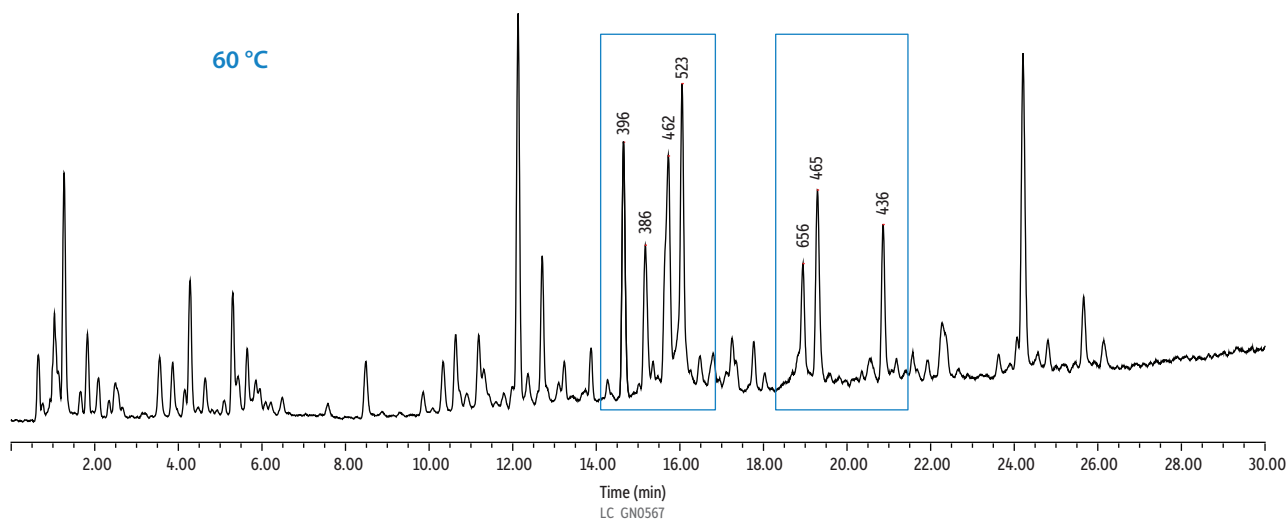


# Peptides on Raptor ARC-18 Poster Figure 4: Effects of Temp.

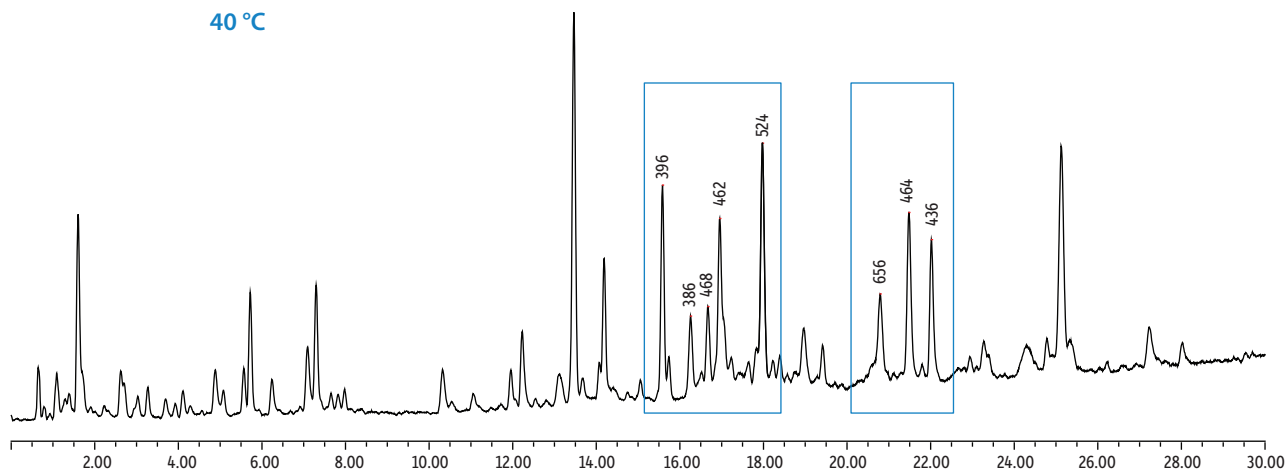

**Column** Raptor ARC-18 (cat.# 9314A12)  
**Dimensions:** 100 mm x 2.1 mm ID  
**Particle Size:** 2.7 µm  
**Pore Size:** 90 Å  
**Sample** Bovine serum albumin (BSA) tryptic digest standard  
**Diluent:** Water + 0.1% formic acid  
**Conc.:** 5 pmol  
**Inj. Vol.:** 2 µL  
**Mobile Phase**  
**A:** Water + 0.1% formic acid  
**B:** Acetonitrile + 0.1% formic acid

Time (min)	Flow (mL/min)	%A	%B
0.00	0.3	98	2
40.0	0.3	75	25
40.1	0.3	98	2
50.0	0.3	98	2

**Detector** MS  
**Ion Mode:** ESI+  
**Mode:** MS scan  
**Instrument** UHPLC