


# Amitriptyline on Ultra IBD, pH 3

## Peak List

1. amitriptyline

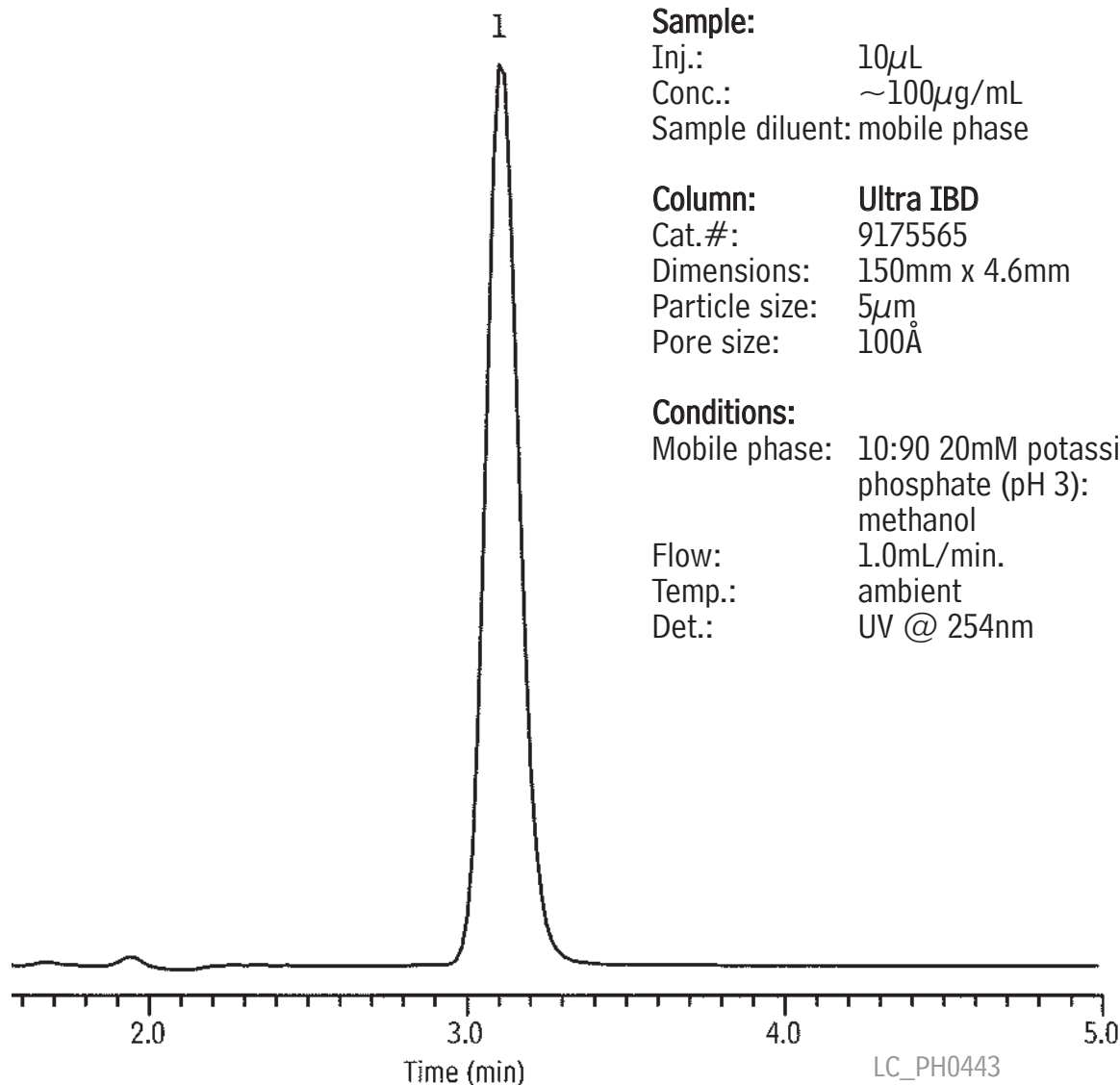

## Sample:

Inj.: 10 $\mu$ L  
Conc.: ~100 $\mu$ g/mL  
Sample diluent: mobile phase

**Column:** Ultra IBD  
Cat.#: 9175565  
Dimensions: 150mm x 4.6mm  
Particle size: 5 $\mu$ m  
Pore size: 100Å

## Conditions:

Mobile phase: 10:90 20mM potassium phosphate (pH 3): methanol  
Flow: 1.0mL/min.  
Temp.: ambient  
Det.: UV @ 254nm