

Amitriptyline on Ultra IBD, pH 7

Peak List

1. amitriptyline


Sample:

Inj.: 10 μ L
Conc.: ~100 μ g/mL
Sample diluent: mobile phase

Column: Ultra IBD
Cat.#: 9175565
Dimensions: 150mm x 4.6mm
Particle size: 5 μ m
Pore size: 100Å

Conditions:

Mobile phase: 20mM potassium phosphate
(pH 7):methanol, 10:90
Flow: 1.0mL/min.
Temp.: ambient
Det.: UV @ 254nm

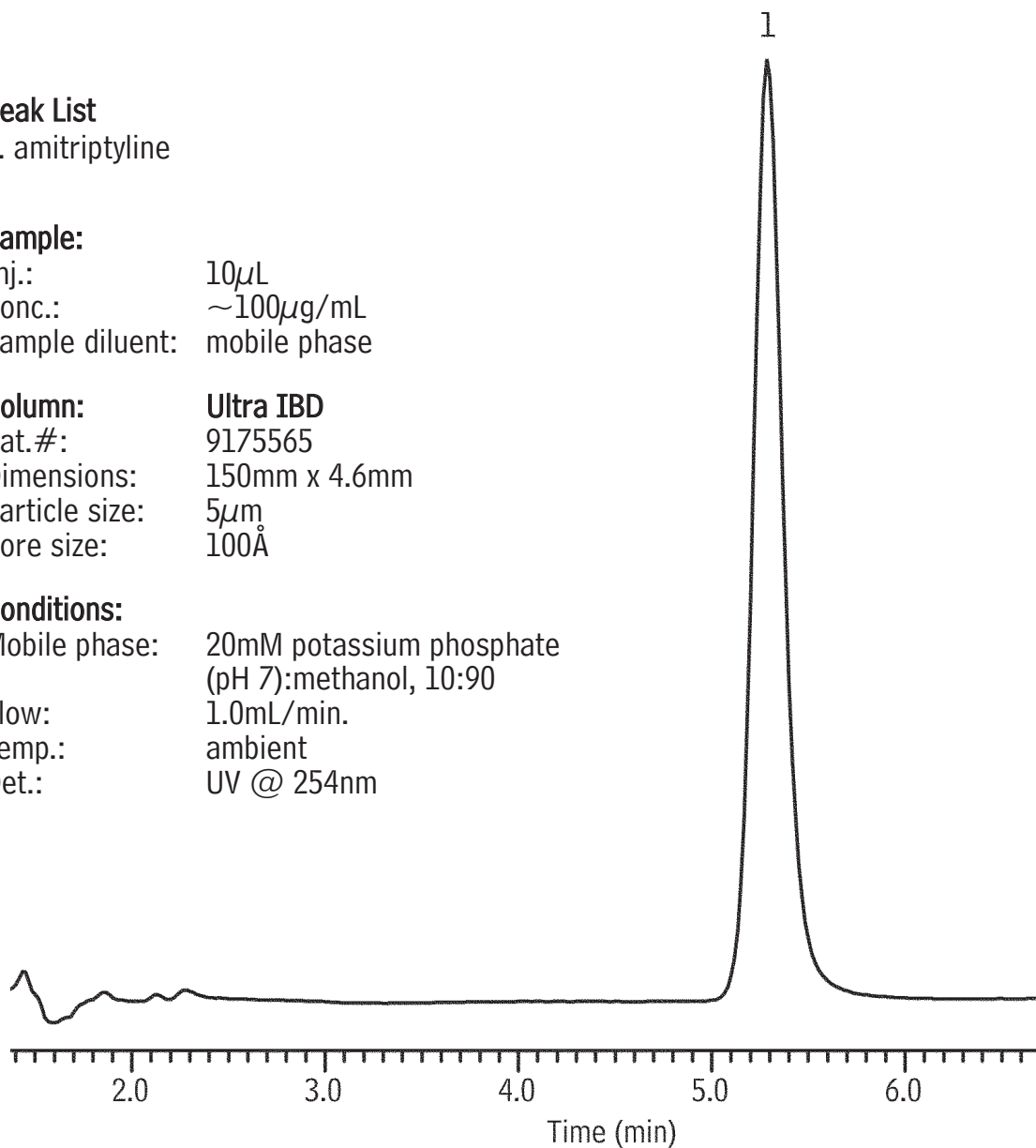

LC_PH0444